

MDR1 C3435T ToolSet™ for LightCycler™ (Pharmakogenetics, P-Glycoprotein)

Lyophilized ToolSet for PCR using the LightCycler™ Instrument. Licensed by Roche Diagnostics GmbH

Order#: MDR1 3435 - 16

1 ToolSet for 16 reactions

Store at 4°C, protected from light.
Exposure to light may especially damage
the Oligotool™ tube (vial with red cap).

For use with LightCycler-DNA Master Hybridization Probes, 10 x conc. (Roche Cat.No.: 2 015 102)

1. ToolSet contents

Vial	Label	Content	Quantity
			MDR1 3435 - 16
1, Red cap	OligoTool	- lyophilized oligos for PCR - contains mutation detection and anchor probe, primers	For 16 tests Dissolved: 50 µL
2, Green cap	Control	- lyophilized Heterozygous DNA	Dissolved: 20 µL
3, Blue cap	Solvent	- to dissolve OligoTool / Control	1000 µL of Solvent

Additional equipment and reagents required but not supplied :
LightCycler-DNA Master Hybridization Probes, 10 x conc. Cat.No.: 2 015 102, including 25mM MgCl₂; LightCycler
instrument, LightCycler capillaries, DNA extraction materials

2. Introduction

2.1. Product overview

ToolSet description This ToolSet is specifically designed for genotyping the C3435T mutation in the MDR1 gene by LightCycler PCR with Melting Curve Analysis. The primer pair and fluorescent detection and anchor probes have been optimized for specific amplification of a 197 bp segment containing the potentially mutated site and optimal genotype discrimination.

Control material Heterozygous control DNA, lyophilized.

Storage of ToolSet and Solutions Store at +4°C when lyophilized, protected from light.
The unopened lyophilized ToolSet is stable at +4°C for 12 months from date of manufacture if protected from light. When dissolved store at +4°C for a maximum of 4 weeks, or at -20°C for longer periods (months), protected from light.
Avoid freezing and thawing.

3. Preparation for LightCycler PCR

Toolset preparation **Dissolve** the content of the **OligoTool** tube (Red Cap) with **50 µl of Solvent**.
Dissolve the content of the **Control** tube (Green Cap) with **20 µl of Solvent**.

1. Before opening tubes, centrifuge them quickly.
2. Add Solvent into OligoTool tube and Control tube as above.
3. Recap tubes, vortex gently.
4. Before opening tubes, centrifuge them quickly.
5. Proceed to Reaction Mix preparation.

Primers ? You don't have to add primers.
Probes ? You don't have to add probes.

Reaction Mix Preparation For 1 (One) reaction, prepare the Reaction Mix as shown in the following table :

Reagent	µL
OligoTool MDR1 3435 -16, dissolved	2.8
Solvent MDR1 3435 -16	10.0
MgCl ₂ 25 mM	1.2 (final 2.5 mM)
Master Hybridization Probes 10x	2
Total Reaction Mix	16
+ Your DNA or Control MDR1 3435 -16	4
Grand Total	20

Use Master Hybridization Probes 10x and MgCl₂ 25 mM from Roche LightCycler-DNA Master Hybridization Probes, 10 x conc. (Roche Cat.No.: 2 015 102, including 25mM MgCl₂).
For multiple reactions, multiply the indicated volumes appropriately.

Positive Control Always run a positive control with the samples. Use the dissolved MDR1 C3435T Heterozygous Control DNA (Green Cap).

Negative control Always run a negative control with the samples. To prepare a negative control, replace the template DNA with Solvent (Blue Cap).

Extraction of genomic DNA You can use different Kits for DNA isolation, either with a manual method or with an automated system. The elution buffers should be salt-free. Example : Roche High Pure PCR Template Preparation Kit (Cat.No. 1 796 828).

Application The **MDR1 gene** encodes the integral membrane protein **P-glycoprotein**, whose function is the energy-driven export of substances from the inside of cells to the outside, including many drugs. The **C3435T** polymorphism correlates with intestinal MDR1 expression and activity. The **T allele** appears to be associated with markedly **lower MDR1** expression compared with the C allele.

Note : This ToolSet was developed for use in life science research only.

4. LightCycler Settings and Experimental Protocol

Denaturation

Cycle Program Data	Value
Cycles	1
Analysis Mode	None
Temperature Targets	Segment 1
Target Temperature (°C)	95
Incubation time (s)	30
Temperature Transition Rate (°/s)	20
Secondary Target Temperature (°C)	0
Step Size (°C)	0
Step Delay (Cycles)	0
Acquisition Mode	None

Amplification

Cycle Program Data	Value		
Cycles	45		
Analysis Mode	None		
Temperature Targets	Segment 1	Segment 2	Segment 3
Target Temperature (°C)	95	55	72
Incubation time (s)	5	5	8
Temperature Transition Rate (°/s)	20	20	5
Secondary Target Temperature (°C)	0	0	0
Step Size (°C)	0	0	0
Step Delay (Cycles)	0	0	0
Acquisition Mode	None	Single	None

Melting Curve Analysis

Cycle Program Data	Value		
Cycles	1		
Analysis Mode	Melting Curves		
Temperature Targets	Segment 1	Segment 2	Segment 3
Target Temperature (°C)	95	40	80
Incubation time (s)	30	120	0
Temperature Transition Rate (°/s)	20	20	0.1
Secondary Target Temperature (°C)	0	0	0
Step Size (°C)	0	0	0
Step Delay (Cycles)	0	0	0
Acquisition Mode	None	None	Continuous

Cooling

Cycle Program Data	Value
Cycles	1
Analysis Mode	None
Temperature Targets	Segment 1
Target Temperature (°C)	40
Incubation time (s)	30
Temperature Transition Rate (°/s)	20
Secondary Target Temperature (°C)	0
Step Size (°C)	0
Step Delay (Cycles)	0
Acquisition Mode	None

LC Program Version and Fluorescence Display Mode

Developed with LC Program Version 3.5 and automatic gain control.
For fluorescence display use F2/F1, or F2 with colour compensation.

5. Typical results

Introduction

Use the Melting Curve program to genotype the human genomic DNA research samples. The melting peaks allow discrimination between the possible genotypes at the **C3435T** mutation site in the **MDR1** gene. Figure 1 shows a typical result obtained with the **MDR1 C3435T ToolSet™** for LightCycler™ :

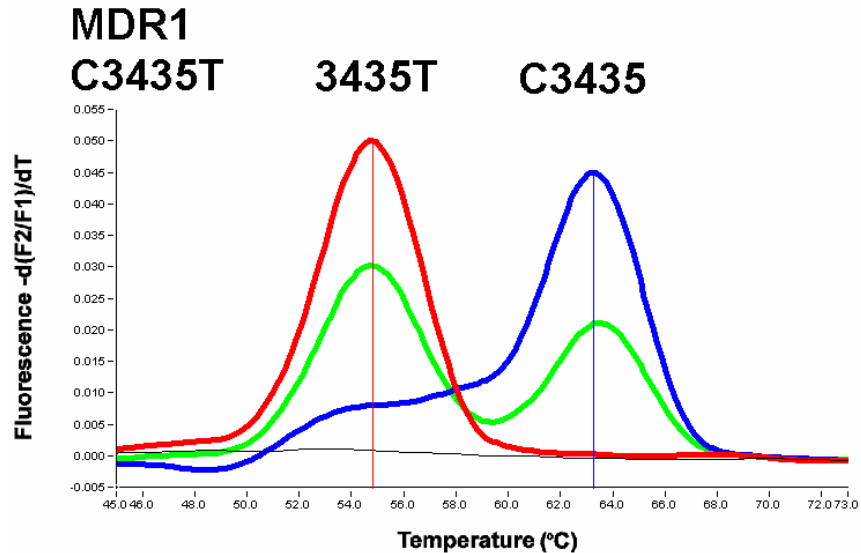


Figure 1 : Melting curve analysis of genotypes at position 3435 of the human MDR1 gene

BLUE : Homozygote **CC 3435** (Wild Type)
GREEN : Heterozygote **C3435T** (HET Control DNA contained in the ToolSet)
RED : Homozygote **3435TT**
BLACK : No DNA Control. **Blue Cursor :** $T_m = 63.3\text{ }^{\circ}\text{C}$; **Red Cursor :** $T_m = 54.8\text{ }^{\circ}\text{C}$

Conditions : LC program version 3.5 with automatic gain setting, No Color compensation, Digital Filter enabled, Degrees to average : 10.0. Calculation Method : Polynomial.

Note : The values for the respective melting temperatures may vary for +/- 2.5 °C between different experiments. The Delta T between the melting peaks for different genotypes may vary +/- 1.0 °C. The MDR1 C3435T ToolSet™ has been developed for and validated with the LightCycler™ and its original accessory materials and reagents. Performance of the Tool-Set with other instruments, accessories and reagents has not been validated by ratiogen.

7. Notices to Purchaser : Licenses and Trademarks, Prohibition of Resale

ratiogen's Genes-4U ToolSets™ for LightCycler™ are produced and marketed under license from Roche Diagnostics GmbH. The purchase of this product does not convey any right for its use in clinical diagnostics. No license to use the PCR process is conveyed expressly or by implication by purchase of this product. The LightCycler technology is licensed from Idaho Technology Inc., Idaho Falls, ID, USA. LightCycler and High Pure are trademarks of members of the Roche group. Resale of ratiogen's Genes-4U products is expressly prohibited.

How to contact ratiogen

E-mail info@ratiogen.com

Internet <http://www.ratiogen.com>